

// **PRODUCT PORTFOLIO**


Veterinary segment



VET LIFT

Veterinary surgery and treatment table

- High-quality scissor-lift table with versatile use in veterinary medicine
- Frame with maximum stability; powder-coated in RAL 7035
- Electrically, battery-powered or hydraulically height-adjustable via foot switch within a wide range from 31–101 cm
- Top (stainless steel or synthetic) with edge and drain (130 × 60 cm), 14 tie-down hooks
- Either stationary with 2 one-side integrated castors (allowing easy movement of the table by lifting on the opposite side) or mobile with 4 lockable double castors
- Table top optionally fixed or tilting (continuously adjustable tilt angle of the tabletop up to 11°) with comfortable one-hand operation, supported by gas spring
- Variant with integrated scale optionally available

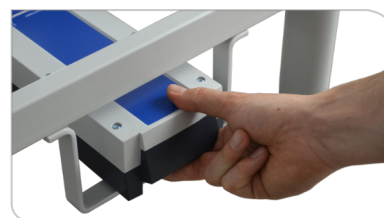

Model variant

Item code (stationary)

Item code. (movable)

VET LIFT, electric, with fixed stainless-steel top	759605	759601
VET LIFT, electric, with tilting stainless-steel top	759606	759602
VET LIFT, electric, with fixed synthetic top	759615	759611
VET LIFT, electric, with tilting synthetic top	759616	759612
VET LIFT, hydraulic, with fixed stainless-steel top	769605	769601
VET LIFT, hydraulic, with tilting stainless-steel top	769606	769602
VET LIFT, hydraulic, with fixed synthetic top	769615	769611
VET LIFT, hydraulic, with tilting synthetic top	769616	769612
VET LIFT, battery-powered, with fixed stainless-steel top	779605	779601
VET LIFT, battery-powered, with tilting stainless-steel top	779606	779602
VET LIFT, battery-powered, with fixed synthetic top	779615	779611
VET LIFT, battery-powered, with tilting synthetic top	779616	779612

Version with integrated scale

Model variant

Item code



VET LIFT, stationary, electric, fixed stainless-steel top, with integrated scale	759607
VET LIFT, stationary, electric, fixed synthetic top, with integrated scale	759617
VET LIFT, stationary, battery-powered, fixed stainless-steel top, w. intg. scale	779607
VET LIFT, stationary, battery-powered, fixed synthetic top, w. integrated scale	779617
VET LIFT, movable, electric, fixed stainless-steel top, with integrated scale	759603
VET LIFT, movable, electric, fixed synthetic top, with integrated scale	759613
VET LIFT, movable, battery-powered, fixed stainless-steel top, w. intg. scale	779603
VET LIFT, movable, battery-powered, fixed synthetic top, w. integrated scale	779613



VET LIFT STAINLESS

Veterinary surgery and treatment table

- Manufactured from high-quality stainless steel
- Electro-chemically polished, smooth, glossy surface (high corrosion resistance, optimal cleaning behavior, significantly reduced tendency to deposits and chemically passive)
- Electrically or battery-powered height adjustment via foot switch from 31–101 cm
- Top (stainless steel or synthetic) with edge and drain (130 × 60 cm), 14 tie-down hooks
- Either stationary with 2 single-side integrated castors (allowing easy movement of the table by lifting on the opposite side) or mobile with 4 lockable twin castors
- Table top optionally fixed or tilting (continuously adjustable tilt angle of the worktop up to 11°) with comfortable one-hand operation, supported by gas spring
- Variant with integrated scale optionally available

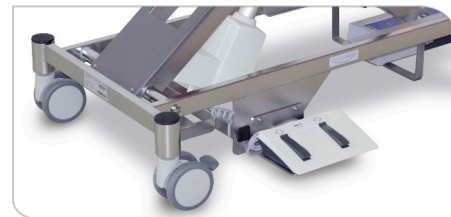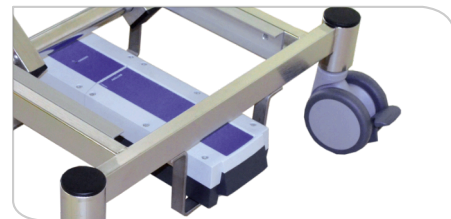

Model variant	Item code
VET LIFT STAINLESS, stationary, electric, with fixed stainless-steel top	659605
VET LIFT STAINLESS, stationary, electric, with tilting stainless-steel top	659606
VET LIFT STAINLESS, stationary, electric, with fixed synthetic top	659615
VET LIFT STAINLESS, stationary, electric, with tilting synthetic top	659616
VET LIFT STAINLESS, stationary, battery-powered, w. fixed stainless-steel top	679605
VET LIFT STAINLESS, stationary, battery-powered, w. tilting stainl.-steel top	679606
VET LIFT STAINLESS, stationary, battery-powered, with fixed synthetic top	679615
VET LIFT STAINLESS, stationary, battery-powered, with tilting synthetic top	679616
VET LIFT STAINLESS, movable, electric, with fixed stainless-steel top	659601
VET LIFT STAINLESS, movable, electric, with tilting stainless-steel top	659602
VET LIFT STAINLESS, movable, electric, with fixed synthetic top	659611
VET LIFT STAINLESS, movable, electric, with tilting synthetic top	659612
VET LIFT STAINLESS, movable, battery-powered, w. fixed stainless-steel top	679601
VET LIFT STAINLESS, movable, battery-powered, w. tilting stainl.-steel top	679602
VET LIFT STAINLESS, movable, battery-powered, with fixed synthetic top	679611
VET LIFT STAINLESS, movable, battery-powered, with tilting synthetic top	679612

Version with integrated scale

Model variant	Item code
VET LIFT STAINLESS, stationary, electric, with fixed stainless-steel top, with integrated scale	659607
VET LIFT STAINLESS, stationary, electric, with fixed synthetic top, with integrated scale	659617
VET LIFT STAINLESS, stationary, battery-powered, with fixed stainless-steel top, with integ. scale	679607
VET LIFT STAINLESS, stationary, battery-powered, with fixed synthetic top, with integrated scale	679617
VET LIFT STAINLESS, moveable, electric, with fixed stainless-steel top, with integrated scale	659603
VET LIFT STAINLESS, moveable, electric, with fixed synthetic top, with integrated scale	659613
VET LIFT STAINLESS, moveable, battery-powered, with fixed stainless-steel top, with integ. scale	679603
VET LIFT STAINLESS, moveable, battery-powered, with fixed synthetic top, with integrated scale	679613

TREND

Veterinary surgery and treatment table

- Made of high-quality steel and aluminum
- Powder-coated surface
- Column-lift system and chassis with maximum stability
- Electrically height-adjustable (230 V) via foot switch from 72–102 cm
- Foot switch operable from both sides
- Choice of stainless-steel or synthetic top (130 × 60 cm) with edge and drain
- One-side integrated castors for easy cleaning
- Continuously adjustable tilt angle of the worktop up to almost vertical
- Comfortable one-hand operation, supported by gas spring
- Special shape of the base frame simplifies power supply via floor socket or energy column
- Chassis housing available in various colors
- Option with integrated heating also available



Version with standard top

Model variant

Item code

TREND column-lift table, electric, with tilting stainless-steel top
TREND column-lift table, electric, with tilting synthetic top

754601cccc
754611cccc

TREND column-lift table, electric, with tilting stainless-steel top (heatable)

754601ccccSH

Version with V-Top

- Removable stainless tops (Total size 143x50 cm)
- Continuously adjustable V-Top (up to 90°) with comfortable one-hand operation supported by gas pressure spring
- Drain gutter and drip tray (removable) as well as longitudinal OP side rails as well as 4 tie-down cleats included



Model variant

Item code

TREND column-lift table, electric, with V-Top stainless-steel top

754531cccc

Color palette (*additional special colors for the powder-coating may be available on request and at extra charge)

Nr. 4002

red-violet

Nr. 1007

daffodil yellow

Nr. 6027

light green

Nr. 5005

signal blue

Nr. 7035

light grey

ECO LIFT

Veterinary surgery and treatment table

- Column lift table with powder-coated surface and chassis with maximum stability
- Hydraulically or electrically height-adjustable (230 V) via foot switch from 72–102 cm
- Foot switch operable from both sides
- Anti-slip adjustable rubber feet for solid stand and compensating potentially uneven floors
- Continuously adjustable tilt angle of the worktop up to almost vertical
- Chassis housing available in various colors
- Option with integrated heating or washing table also available



Version with standard top

Model variant	Item code
ECO LIFT column-lift table, electric, with tilting stainless-steel top	651600cccc
ECO LIFT column-lift table, electric, with tilting synthetic top	651610cccc
ECO LIFT column-lift table, hydraulic, with tilting stainless-steel top	661600cccc
ECO LIFT column-lift table, hydraulic, with tilting synthetic top	661610cccc
ECO LIFT column-lift table, electric, with tilting stainless-steel top (heating)	651600ccccSH
ECO LIFT column-lift table, hydraulic, with tilting stainless-steel top (heating)	661600ccccSH
ECO LIFT column-lift table, electric, with washing table	651640cccc
ECO LIFT column-lift table, hydraulic, with washing table	661640cccc

Version with V-Top

- Removable stainless tops (total size 143x50 cm)
- Manually adjustable tilting of tabletop (up to approx. 45°) as well as adjustable V-Top (up to 90°)
- Drain gutter and drip tray (removable) as well as longitudinal OP side rails included



Model variant	Item code
ECO LIFT column-lift table, electric, with stainless-steel V-Top	651530cccc
ECO LIFT column-lift table, hydraulic, with stainless-steel V-Top	661530cccc
ECO LIFT column-lift table, electric, with stainless-steel V-Top (heating)	651530ccccSH
ECO LIFT column-lift table, hydraulic, with stainless-steel V-Top (heating)	661530ccccSH

Color palette (*additional special colors for the powder-coating may be available on request and at extra charge)

Nr. 4002 red-violet	Nr. 1007 daffodil yellow	Nr. 6027 light green	Nr. 5005 signal blue	Nr. 7035 light grey
-------------------------------	------------------------------------	--------------------------------	--------------------------------	-------------------------------

BASE LIFT

Veterinary surgery and treatment table

- Robust veterinary surgery and treatment table made of high-quality materials and components
- Either premium stainless-steel or synthetic top with edge and drain (130 × 60 cm)
- 14 tie-down hooks included
- Height-adjustable from 54–100 cm, either electrically (power-cord) or battery-powered
- Table top optionally fixed or tilting via gas spring (continuously adjustable up to approx. 30°) for comfortable one-hand operation
- Foot switch for height adjustment
- Either stationary on anti-slip feet or movable on easy-to-maneuver double castors, each individually lockable
- Extensive range of accessories available on request
- Elegantly pressed table top; full surface underneath is supported by wooden core



Model variant

Item code

BASE LIFT, stationary, electric, with fixed stainless-steel top	638605
BASE LIFT, stationary, electric, with tilting stainless-steel top	638600
BASE LIFT, stationary, electric, with fixed synthetic top	638615
BASE LIFT, stationary, electric, with tilting synthetic top	638610
BASE LIFT, stationary, battery-powered, with fixed stainless-steel top	648605
BASE LIFT, stationary, battery-powered, with tilting stainless-steel top	648600
BASE LIFT, stationary, battery-powered, with fixed synthetic top	648615
BASE LIFT, stationary, battery-powered, with tilting synthetic top	648610
BASE LIFT, moveable, electric, with fixed stainless-steel top	638606
BASE LIFT, moveable, electric, with tilting stainless-steel top	638601
BASE LIFT, moveable, electric, with fixed synthetic top	638616
BASE LIFT, moveable, electric, with tilting synthetic top	638611
BASE LIFT, moveable, battery-powered, with fixed stainless-steel top	648606
BASE LIFT, moveable, battery-powered, with tilting stainless-steel top	648601
BASE LIFT, moveable, battery-powered, with fixed synthetic top	648616
BASE LIFT, moveable, battery-powered, with tilting synthetic top	648611

ECO VET

Veterinary surgery and treatment table

- Electrically (power-cord) or battery-powered height adjustment from 54–100 cm
- Either premium stainless-steel or synthetic top with edge and drain (130 × 60 cm)
- 14 tie-down hooks included
- Table top optionally fixed or tilting via mechanical adjustment (continuously adjustable up to approx. 30°)
- Foot switch for height adjustment
- Either stationary on anti-slip feet or movable on easy-to-maneuver, individually lockable castors
- Extensive range of accessories available on request
- Elegantly pressed table top; full surface underneath is supported by wooden core



Model variant	Item code
ECO VET lift table, stationary, electric, with tilting stainless-steel top	758600
ECO VET lift table, stationary, electric, with fixed stainless-steel top	758605
ECO VET lift table, stationary, electric, with tilting synthetic top	758610
ECO VET lift table, stationary, electric, with fixed synthetic top	758615
ECO VET lift table, stationary, battery-powered, with tilting stainless-steel top	778600
ECO VET lift table, stationary, battery-powered, with fixed stainless-steel top	778605
ECO VET lift table, stationary, battery-powered, with tilting synthetic top	778610
ECO VET lift table, stationary, battery-powered, with fixed synthetic top	778615
ECO VET lift table, moveable, electric, with tilting stainless-steel top	758601
ECO VET lift table, moveable, electric, with fixed stainless-steel top	758606
ECO VET lift table, moveable, electric, with tilting synthetic top	758611
ECO VET lift table, moveable, electric, with fixed synthetic top	758616
ECO VET lift table, moveable, battery-powered, with tilting stainless-steel top	778610
ECO VET lift table, moveable, battery-powered, with fixed stainless-steel top	778606
ECO VET lift table, moveable, battery-powered, with tilting synthetic top	778611
ECO VET lift table, moveable, battery-powered, with fixed synthetic top	778616

PURE

Veterinary surgery and treatment table

- Chassis manufactured from highest-grade steel with maximum stability; powder-coated RAL 9007
- Particularly hygienic housing made of fully enclosed, high-quality stainless steel with protective sealing lips
- Electrically height-adjustable via foot switch
- Minimum height of 11 cm or 18 cm (including stand foot)
- Maximum height of 112 cm or 119 cm (including stand foot)
- Top either stainless steel or synthetic; optionally with integrated scale
- Either free-standing with stand foot or wall-mounted for maximum stability
- Table top optionally fixed or tilting



Version with stand foot

Model variant	Item code
PURE treatm. & surgery table, electric, with fixed stainless-steel top	655600
PURE treatm. & surgery table, electric, with tilting stainless-steel top	655602
PURE treatm. & surgery table, electric, with integrated scale & stainless top	655603
PURE treatm. & surgery table, electric, with fixed synthetic top	655610
PURE treatm. & surgery table, electric, with tilting synthetic top	655612
PURE treatm. & surgery table, electric, with integrated scale & synthetic top	655613

Version with wall-mounting

Model variant	Item code
PURE treatm. & surgery table, electric, with fixed stainless-steel top	655600-WM
PURE treatm. & surgery table, electric, with tilting stainless-steel top	655602-WM
PURE treatm. & surgery table, electric, with integrated scale / stainless top	655603-WM
PURE treatm. & surgery table, electric, with fixed synthetic top	655610-WM
PURE treatm. & surgery table, electric, with tilting synthetic top	655612-WM
PURE treatm. & surgery table, electric, with integrated scale / synthetic top	655613-WM

OP-System

X-ray compatible special OR table for C-arm

- High-end veterinary special OR table for ophthalmology and advanced orthopedics
- Designed for complex procedures
- Frame of maximum stability, powder-coated in RAL 7035
- Surgery also possible in a seated position
- Height-adjustable from 59–99 cm, either electric or battery-powered, via hand control (optional joystick or multi-function foot switch)
- Movable on centrally and directionally lockable, electrically conductive double castors; activated via pedal for optimal maneuverability
- Electrically adjustable laterally (up to 40°) and longitudinally in both directions (up to 15°)
- Extendable top (X-ray compatible), for example for ophthalmology or C-arm imaging
- Standard scope of delivery: 4 tie-down cleats, OP rails all around (4 rails in total)



Model variant

Item code

OP-System, moveable, electric, with sliding top (X-ray compatible)
OP-System, moveable, battery-powered, with sliding top (X-ray compatible)

652541
672541

Version with V-Top*

- High-quality, removable stainless-steel tops (total size 143 × 50 cm)
- Continuously adjustable V-Top (up to 90°) with comfortable one-hand operation, supported by gas spring
- Extendable drain gutter, 4 tie-down cleats as well as longitudinal OR rails included as standard



**In contrast to the version with extendable top: not X-ray compatible*

Model variant

Item code

OP-System, moveable, electric, with stainless-steel V-Top
OP-System, moveable, battery-powered, with stainless-steel V-Top

652531
672531

OP-System V3 X-Top

Special OR table suitable for laparoscopy

- High-end veterinary special OR table for ophthalmology and advanced orthopedics
- Column-lift chassis with maximum stability
- Removable and re-insertable laparoscopy cut-outs on both sides (individually), allowing closer working at the patient
- Developed for complex and demanding procedures – surgery also possible in seated position
- Electrically height-adjustable from 64–104 cm
- Electrically adjustable laterally (up to 40°)
- V-position up to 90° via comfortable one-hand operation, supported by gas spring
- Longitudinal adjustment electric up to 10° in both directions for (anti-)Trendelenburg
- Extensive space underneath the table top compared with the previous version thanks to optimized design
- Movable on centrally and directionally lockable, electrically conductive double castors; activated via pedal for optimal maneuverability
- Removable, easy-to-clean drain gutter with collecting drip tray for fluids
- Standard scope of delivery: 4 tie-down cleats, OP rails all around (4 rails in total), 2 removable laparoscopy inserts



Model variant

Item code


OP-System V3 X-Top, moveable, electric, with removable laparoscopy cut-outs
OP-System V3 X-Top, moveable, battery-powered, with removable laparoscopy cut-outs

652551
672551

AERON

X-ray compatible OR table for C-arm

- Electrically height-adjustable from 80–110 cm
- Electrically tilting longitudinally 15°,
- Electrically tilting laterally 20°
- X-ray translucent synthetic top (X-ray area depending on chosen variant), extremely resistant
- Suitable for arthroscopy in small animals (extensive accessory range on request)
- With battery or power-cord connection
- Mobile on double castors
- Fixation on individually lockable, anti-slip rubber feet for excellent stability
- Lifting capacity: 120 kg
- Version with V-shaped top features individually removable frame halves: to enlarge the X-ray area (centrally split frame) up to almost entire tabletop (incl. Sides)

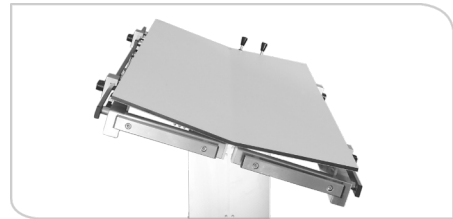

Standard equipment

- OR rails (all around)
- 4 tie-down cleats
- Joystick control

Model variant

Item code

AERON, battery-powered, with V-shaped synthetic top (X-ray compatible)
AERON, electric, with V-shaped synthetic top (X-ray compatible)

776002
756002

AERON, battery-powered, with flat synthetic top (X-ray compatible)
AERON, electric, with flat synthetic top (X-ray compatible)

776522
756522

Special Equipment

Ultrasound table

- Ideal for complex, cardio-logical ultrasound examinations with excellent patient comfort
- Water-repellent, anti-static upholstery made of hard-wearing synthetic leather
- Optimal cleaning behavior and greatly reduced tendency to deposits



Model variant


Ultrasound table (powder-coated) with triangular cut-out
Ultrasound table (stainless steel) with triangular cut-out and insert

Item code

602.640.0-PB
602.640.1

Overlay mat (seamless)

- Seamless, hygienic and antibacterial special coating
- Extremely durable and highly scratch-resistant and waterproof, as well as bite-resistant
- Developed in cooperation with leading experts of the University of Veterinary Medicine, Vienna

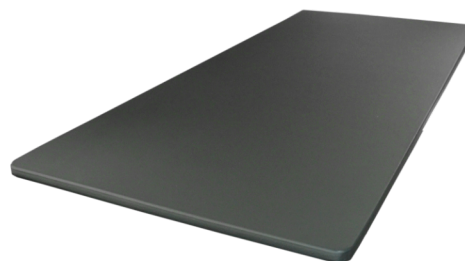

Model variant

Overlay mat black (seamless) 123 x 53 cm (for tops 130 x 60 cm)
Overlay mat black (seamless) 113 x 53 cm (for tops 120 x 60 cm)

Item code

70.123.053.02-P
70.113.053.02-P

Folding table (stainless steel)

- High-quality, polished stainless steel
- Wall-mounted and fully fold-able; incl. full-surface stainless-steel wall protection
- Flat top (90 x 50 cm) optionally in synthetic or stainless steel (load capacity up to 100 kg)
- OP rail for comfortable mounting of accessories



Model variant

Folding table V2 with flat stainless-steel top
Folding table V2 with flat synthetic top

Item code

603.090.050.S
603.090.050.C

Washing table

- Grid top in 4 removable sections
- Mounting for tap fitting and drain
- Dimensions: 146 x 65 x 85 cm
- Made of high-quality stainless steel with 10 cm deep basin (special sizes only on request)



Model variant

Washing table stainless steel (stationary)
ECO LIFT washing table stainless steel, electrically height-adjustable
ECO LIFT washing table stainless steel, hydraulically height-adjustable

Item code

71-136.090-W
651.640cccc
661.640cccc

Further items can be found on our website and may additionally be available on request.

Accessories & Spare parts

Steel treatment table (stainless steel)

- Stationary examination and treatment table made of high-quality stainless steel
- Choice of stainless-steel or synthetic top, seamless with edge and drain
- Incl. 14 tie-down hooks
- Stable stand due to robust stainless-steel legs



Model variant	Item code
Stainless steel treatment table (stationary) with stainless-steel top	71-136.090
Stainless steel treatment table (stationary) with synthetic top	71-136.090.K

Cleaning agent for synthetic tops

- Anti-static cleaner and care product for all synthetic tops of our tables
- Prevents electrostatic charging and dust
- Ready-to-use solution, avoids stress cracking and corrosion damage



Model variant	Item code
Cleaner (antistatic) for synthetic tops 250 ml	Z.KSR-AS-0250

Accessoires

Model variant	Item code
OP rails all around (set of 4 pcs.)	602600
OP rails (piece-wise, or pair-wise)	602596, 602597, 602598, 602599
Tie-down cleat	602601
Infusion pole	602602
Adjustable side supporters (1 pair)	602604
Narcotic bow	602605
Leg support with fixation clamp	602606
Tilting steel dish with grid, incl. insert plate	602607
Drawer for X-ray cassette*	602608
OP-rail extension 200 mm (1 pair)	602610
Instrument tray 300 x 200 mm	602611
Arm rest for ophthalmology (1 pair)	602614
Top extension 50 x 30 cm	602615
Top extension 30 x 30 cm	602616
Insert for OP-System V2 (new)*	602.635-OPSN
Lateral positioners with 3 different types of cushions, or with toothed joint	602641, 602646
Top extension 90 x 30 cm, tiltable 90°, X-ray translucent	612232

* Exclusively compatible with OP-System.

Otherwise, accessories are generally compatible only with our OR rails / supply rails.

Spare parts

Model variant	Item code
Spare battery	610000
Charger for battery	610001
Foot switch	880000
Joystick for OP-System	605000
Castors set (d=75mm) for ECO LIFT or ECO VET	892001
Hand switch 3-channel for OP-System	ZHB83
Power cord with EU-plug	ZNEKA-EU-G-G

Further spare parts available on request. For any additional spare parts required for our veterinary portfolio, please contact us – ideally already stating the serial number of the device concerned.

www.panno-med.at

QUALITY & DESIGN - MADE IN AUSTRIA

PANNOMED