



Excellence in illuminating clinical care begins with precision at the millimeter level

An Ideal Surgical Light:

- **Accurately reproduces the natural colors and subtle differences in tissue**, enabling veterinarians to make precise clinical decisions and perform safe, meticulous procedures.
- **Provides homogeneity and stable illumination**, reduces the need for frequent light adjustments that may interfere with surgical workflow.
- **Reduces eye fatigue**, promoting surgeons' eye health to assist in focused and composed high-intensity cases.
- **Adjusts brightness transitions in accordance with visual physiology**, reducing discomfort in animals caused by sudden exposure to intense light.



Q | Search Mindray Animal Medical

www.mindrayanimal.com
service@mindrayanimal.com

P/N:ENG-VeLED Series-210285X6P-20250505
©2025 Shenzhen Mindray Animal Medical Technology Co., LTD. All rights reserved.

| Model | Surgical Lights | | Examination Lights | |
|--------------------------------|---|---|--------------------------------------|---------------|
| | VeLED S5 | VeLED S3 | VeLED E5 | VeLED E3 |
| Illumination (Lux) | 160000 | 100000 | 70000 | 60000 |
| Light Field Diameter (mm) | 170-250 | 170-250 | 170 | 170 |
| 60% Depth of Illumination (mm) | 600 | 600 | 1200 | 1200 |
| Light Head Size (mm) | 500 | 500 | 290 | 290 |
| Mounting Options | Ceiling : Single & Dual-head Wall/Mobile | Ceiling : Single & Dual-head Wall/Mobile | Ceiling : Single-head Wall/Mobile | Wall/Mobile |
| Arm Configurations | Spring Arm | Spring Arm | Spring Arm | Gooseneck Arm |
| Brightness Control Handle | Standard | Optional | Optional | Optional |
| Endo Mode | √ | √ | × | × |

For detailed parameter configurations, please consult Mindray Animal Medical's authorized distributor



VeLED Series LED Surgical / Examination Lights

Steel may hold the scalpel
but light holds the soul



CRI/ **99**
Accurate Color Restoration

R9/ **97**
Enhanced Vessel Recognition

Color Temperature / **4500**
Visual Strain Reduction

60% Depth of Illumination / **600**
Clear Visualization of Deep Cavity

Exceptional Illumination Eye-friendly Comfort



High Performance & Low Fatigue

Superior color rendering and illumination performance
Reduces visual fatigue during prolonged surgical procedures

Multi-patch Superposition Technology (MPST)

- Delivers consistently uniform and circular light across the surgical field
- Maintains stable illumination even when part of the light source is obstructed

Electronic Focusing

- Inner and outer ring beams converge at 0.8–1.2 m without dispersion
- Reduces need for frequent light adjustments, allowing surgeons to stay focused on procedures



Adjustable Light Field Diameter

- Light field diameter adjustable from 170 to 250mm
- Designed to suit a wide range of surgical cases - from fine incisions to broad operative fields

Stable Articulation & Freely Suspended Positioning

- Maintains precise positioning without drifting or sagging over long-term use
- Enables fast and stable angle and position adjustments during surgery

Dual-Function Handle

- Optional control handle integrates both positioning and brightness adjustment
- Enhances surgical efficiency and workflow fluidity

Endo Mode

- Provides soft, low-illuminance lighting optimized for endoscopic procedures
- Prevents screen glare and maintains a clear visual field



Single-handed operation Full-room flexibility



Integrated Lamp Head

- Screw-free
- Antimicrobial coating
- Effortless routine care


Compact Agility Powerful Functionality



225mm
Ultra-compact
Flange Tube

Ceiling Height Adaptability


- Engineered to suit ceiling heights as low as 2.4 meters
- Ensures a professional, efficient surgical layout even in compact veterinary settings

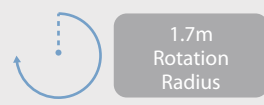

1.7m
Rotation
Radius

Compact Space Solution

- Specifically designed for veterinary clinics with limited space
- Easily fits into surgical areas as small as 5 m², offering flexible and functional installation

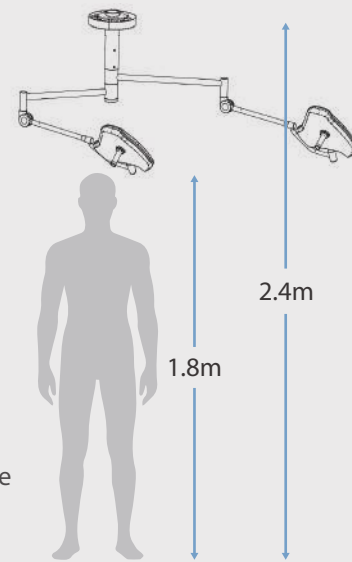

A 5 m² surgical area is one that is specifically configured to support a rotational radius of 1.7 meters

Crafted to Perfection Trusted Quality

- Rigorously tested under extreme conditions to last over 10 years of reliable performance

25
Disinfectants

IP54
The Waterproof and
Dustproof Rating of Lamp Head

10,000+
Aging Tests

20,000
Operation Tests



100,000 Hours
Service of Life

40,000+
Joint Rotation Tests

Low EMI
Minimizes Interference with
Ultrasound Imaging



Multiple mounting options for
flexible and numerous installation options